


Philips HeartStart OnSite AED

Patented Quick Shock feature allows the OnSite to typically deliver a shock within 8 seconds after CPR.¹



1

Pull the handle on Cartridge Cover and press the green On/Off button, which activates voice instruction and visual icons

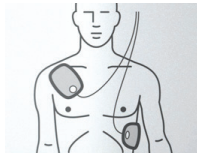


2

Place the pads on the patient as directed.



3

When advised by the device, press the orange Shock button.

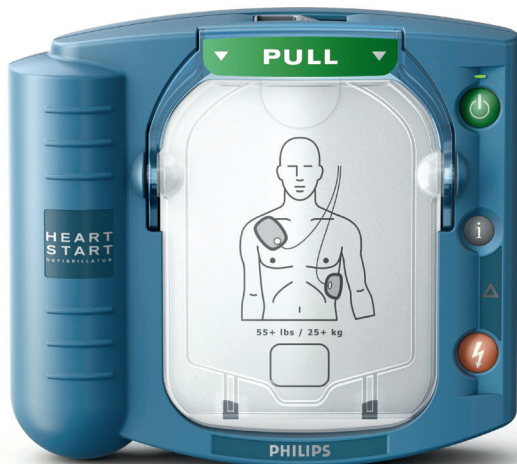

PHILIPS

Automated
External Defibrillator

HeartStart OnSite

Side by side.
Step by step.

Philips HeartStart OnSite AED

1. Nichol, G., Sayre, M. R., Guerra, F., & Poole, J. (2017). Defibrillation for Ventricular Fibrillation: A Shocking Update., 70(12), 1496-1509. *Journal American College of Cardiology* doi:10.1016/j.jacc.2017.07.778

Advanced technology. Proven therapy.

Philips HeartStart OnSite AED

The Philips HeartStart OnSite assists you through the process of treating a victim of suspected sudden cardiac arrest (SCA) and is the only AED available over the counter. The OnSite AED provides practically real-time guidance through step-by-step voice commands and CPR guidance.

- Includes features to help guide the treatment of sudden cardiac arrest with easy setup, clear voice commands and real time metronome
- Arrives virtually ready to use. With the Ready-Pack configuration, the OnSite AED is positioned inside the carry case, with Adult SMART Pads Cartridge and battery already installed and with a spare Adult SMART pads cartridge in place
- Guides you through a cardiac emergency with a simple, step-by-step process, adaptive instructions and intelligent sensors to help deliver therapy
- Use on infants and children under 25 kg or 55 lbs or 0-8 years old, and adults and children over 25 kg or 55 lbs or greater than 8 years old
- Senses when the special Infant/Child SMART Pads cartridge is installed, and automatically adjusts CPR instructions and shock energy
- Can be converted to a trainer with installation of training pads cartridge
- Conducts a series of automatic self-tests daily, weekly and monthly, to check pad readiness and verify functionality and calibration of circuits and systems



Ready to act. Ready to go.

Designed for the ordinary person in the extraordinary moment, Philips HeartStart OnSite AED is ready to act and virtually ready to go. It allows anyone with little or no training to treat the most common cause of sudden cardiac arrest (SCA) by delivering a shock quickly and effectively, wherever SCA happens.

Start quickly. Treat confidently.

With access to the right equipment and support, you can help save a life. The OnSite AED guides you through the process of treating a victim of suspected sudden cardiac arrest. The OnSite AED provides practically real-time guidance through step-by-step voice commands and CPR guidance.



Smart prompts. Simple steps.

We've equipped OnSite with integrated SMART Pads that will provide feedback to the AED so it can adapt its voice instructions to your actions and your pace. The system won't announce the next step until you are ready. Prompts are repeated and rephrased if needed and include additional instruction to aid understanding.