

Baxter

ProBP 2000

Digital Blood Pressure Device,
2000 PLUS-CE Model



Choose our **ProBP 2000 PLUS** device if you're looking to upgrade blood pressure readings from a manual to digital approach.

The ProBP 2000 PLUS device Advantage —

New design comes with a BP averaging feature, to help give a broader view of the patient's true blood pressure than a single reading.



Pre-programmed
three-blood pressure
averaged reading



One button touch
triggers the device
to start capturing



Memory recall
button to review
last reading



SureBP technology offers great advantages in speed and patient comfort...clinicians can have added confidence in this technology¹."

We've Combined Two of Our Technologies to Support a Streamlined Workflow for You and Your Patients.

1 SureBP Blood Pressure Technology


Our **SureBP** inflation-based algorithm is used across **Welch Allyn** vital signs devices, helping clinicians in a range of care settings capture automated, fast and accurate blood pressure readings.

- Blood pressure readings can be captured in about 20 seconds², helping address patient comfort.
- Meets IEC standard for blood pressure accuracy.
- Motion indicator signals to the user when motion is present during the measurement, which may result in an inaccurate reading.

FAST AND AUTOMATED BLOOD PRESSURE READINGS IN A VARIETY OF SET-UPS:



Wall Mount Option with Bracket
(Integrated Wall System sold separately)

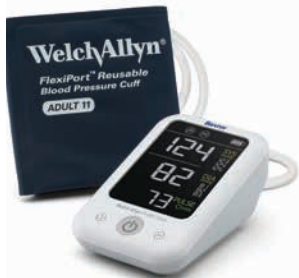

Hand-held

2 FlexiPort Blood Pressure Cuff Technology

Using the correct cuff size is an important factor in obtaining accurate readings. **Welch Allyn FlexiPort** cuffs come in a variety of sizes to help address a wide range of your patient population.

- Our unique 'snap-fit', single point connection design was developed to reduce the need for traditional cuff tubes and connectors.
- Our **ProBP 2000** models are compatible with **FlexiPort** Reusable cuffs, ranging from child up to large adult and thigh (15-55 cm range).



Mobile Stand

ITEM	SPECIFICATION
External Dimensions, Approximate	Width: 94 mm (3.70 in) Length: 142 mm (5.59 in) Height: 66 mm (2.60 in)
Net Weight, Approximate (Excluding Batteries)	Model 2000 PLUS: 354 g (0.78 lbs)
Display Mode	Digital LCD, vertical alignment (VA), 68 mm x 90 mm, pressure and pulse displayed
Power supplies and their AC plug-adaptor modes	Optional
Battery-powered (DC) Mode	6 V, 4 AA alkaline batteries
Measurement model	Oscillometric testing mode
Measurement Range	Rated cuff pressure: 0 mmHg – 300 mmHg (0 kPa – 40 kPa) Systolic Measurement Pressure: 50 – 260 mmHG Diastolic Measurement Pressure: 25 – 220 mmHG Pulse Value: 40 to 199 beats/minute
Accuracy	Pressure: $\pm 0.4\text{kPa}$ (3 mmHg); Pulse Value: $\pm 4\%$
Product Life Expectancy	5 years or 15,000 measurements (whichever comes first)

PARTS AND ACCESSORIES

CONFIGURATIONS

2000 PLUS-CE: ProBP 2000 PLUS Digital Blood Pressure Device with BP averaging and memory recall, Includes Size-11 Adult, Reusable FlexiPort Cuff; 4 x AA Batteries; Quick Reference Guide; power supply UK.

FLEXIPORT BLOOD PRESSURE CUFFS

REUSE-09:	Child (15–21 cm)
REUSE-10:	Small Adult (20–26 cm)
REUSE-11:	Adult (25–34 cm)
REUSE-11L:	Adult Long (25–34 cm)
REUSE-12:	Large Adult (32–43 cm)
REUSE-12L:	Large Adult Long (32–43 cm)
REUSE-13:	Thigh (40–55 cm)

ACCESSORIES

107041:	Power Supply (includes US, EU, UK, AU adapters)
107194:	Basket with holster kit
107201:	Holster Kit/Wall Mount Bracket Accessory
107195:	Mobile Stand Kit

To learn more please contact with your Baxter sales representative.

References

1. Alpert, Bruce S., Clinical evaluation of the Welch Allyn SureBP® algorithm for automated blood pressure measurement. Blood Pressure Monitoring, Vol. 12 No. 4, August 2007.
2. Baxter Data on File

The ProBP 2000 Digital Blood Pressure Device is intended for use in measuring blood pressure and pulse rate in patients at least 3 years of age or older with arm circumferences between 15 cm to 55 cm (approximately 5.9 to 21.7 inches). Medical device (2017/745): Class IIa. Notified Body: TÜV SÜD (CE0123). Manufacturer: Guangdong Transtek Medical Electronics Co., Ltd, Zone A, No. 105 Dongli Road, Torch Development District, Zhongshan, 528437, Guangdong, China. EU Authorised Rep.: MDSS-Medical Device Safety Service GmbH, Schiffgraben 41, 30175 Hannover, Germany. This medical device is a regulated health product which, pursuant to such regulation bears a CE mark. Baxter recommends that you carefully read the detailed instructions for safe and proper use included in the documents accompanying the medical devices. The personnel of healthcare establishments are responsible for the proper use and maintenance of these medical devices.

Hillrom.Baxter.eu/en-GB

Baxter Healthcare Ltd
Spur A, 2nd Floor, Reading International Business Park, Basingstoke Road, Reading, RG2 6DH

Baxter, FlexiPort, ProBP, ProBP 2000, ProBP 2000 PLUS, SureBP and Welch Allyn are trademarks of Baxter International Inc. or its subsidiaries.

UKI-FLC134-260001 (v1.0) 03/2026